


## **For operation of ALEC software by Windows Vista**

- 1) Verify whether UAC becomes effective or not.

The control function of user account is reinforced, and so, in case the application software is installed under the off condition of User Account Control(UAC), the installation may be failed because the records of ActiveX, etc. is not normally executed. Due to the above reason, it is necessary to surely turn UAC on at the time of installation of application software.

### **How to make UAC be on**

1. Select “Start” - “Control Panel” – “User Account” in order.
2. Click “Turn User Account Control on or off”.

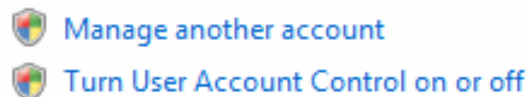

3. Check a mark on “Unse User Account Control(UAC) to help protect your computer”.




3. Re-launch the computer after clicking “OK”.

The installation of ALEC software is able to be normally done after turning UAC on by the above way.


A security mark is added on short-cut of software installed.

2) Setup to compatibility mode of RS232-C EXE


Launch **Explore** and open “Install Folder”. And of the install folder, right-click on the “RS232C EXE”

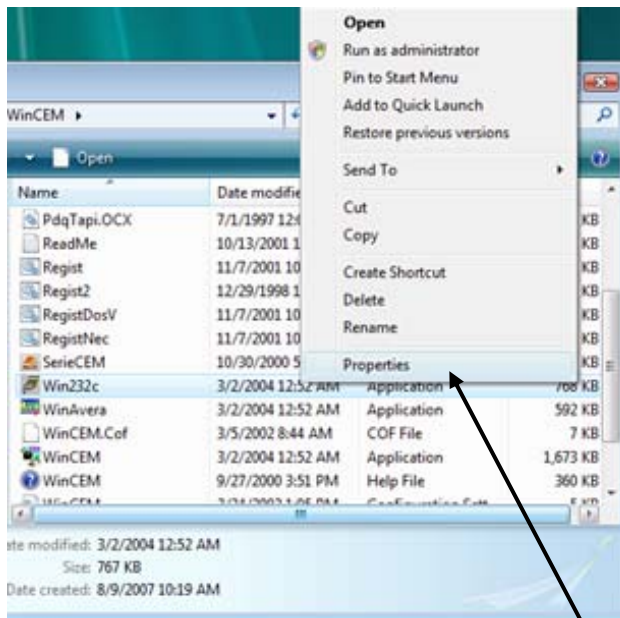

Do drop-down by right-clicking, and select “Properties”.

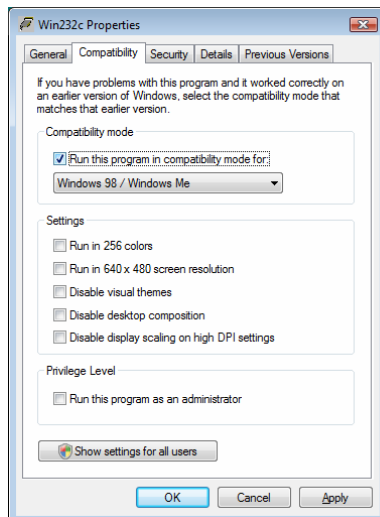


Click a tub of Compatibility and carry on the setup in accordance with the following indication.

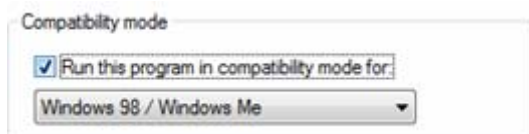

Press “OK” after setup.

### 3) Registry registration

Open “Command Prompto” from “Accessories” and select “Run an administrator”

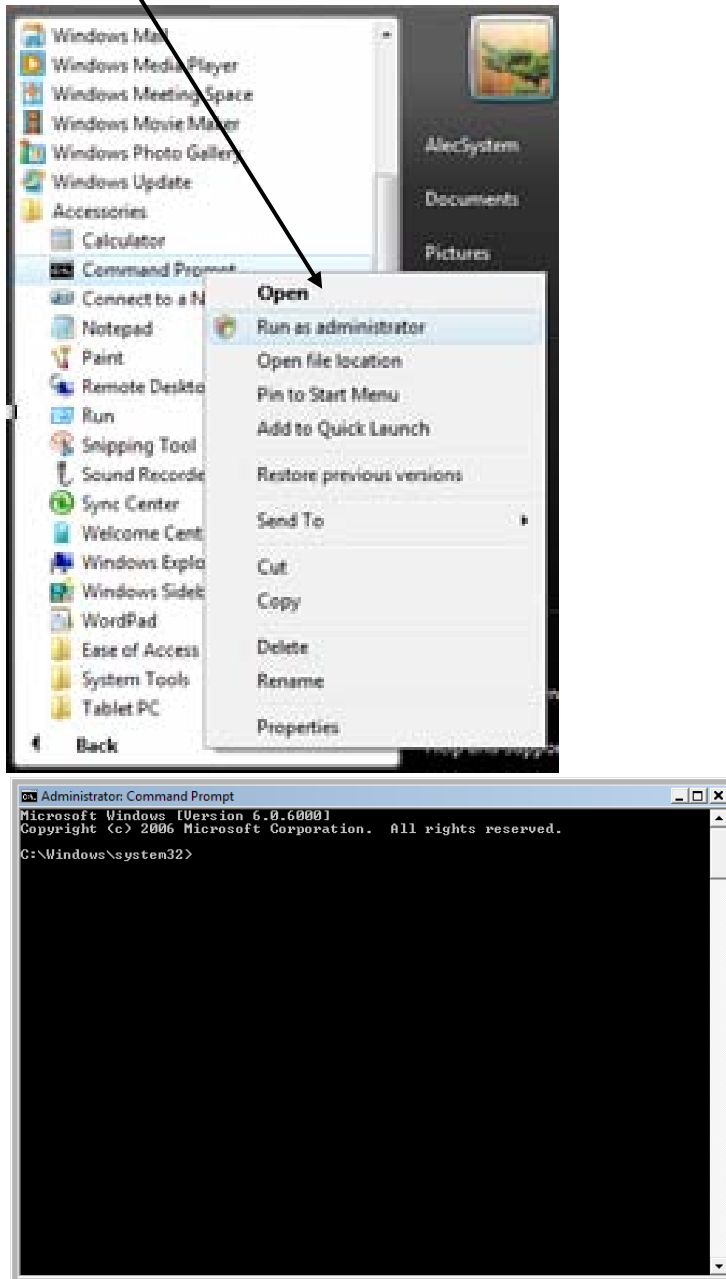

On the “Command Prompto”

CD¥

CD Win XX (XX represents “Promgram name to be used”)

Regist2 C:¥windows

When “Task Completed” is appeared, “Registry” has been completed.